

Newborn bed BeLLY


Newborn bed Belly is an ideal solution that provides the closest contact between a mother and a child. The bed design offers a wide range of positioning, including height adjustment, which enables the bed to be placed over all common beds. The mother can thus care of her child directly from her bed.

Braked wheels, very good stability in all positions, removable bath and its protective frame make the bed perfectly safe. Thanks to the transparent material of the bath you permanently know what the child is doing. A label with child's name can be simply inserted into a transparent pocket on the front of the bed.

Newborn bed Belly is designed for children up to 4 months of age.



Practical
Variable
Safe

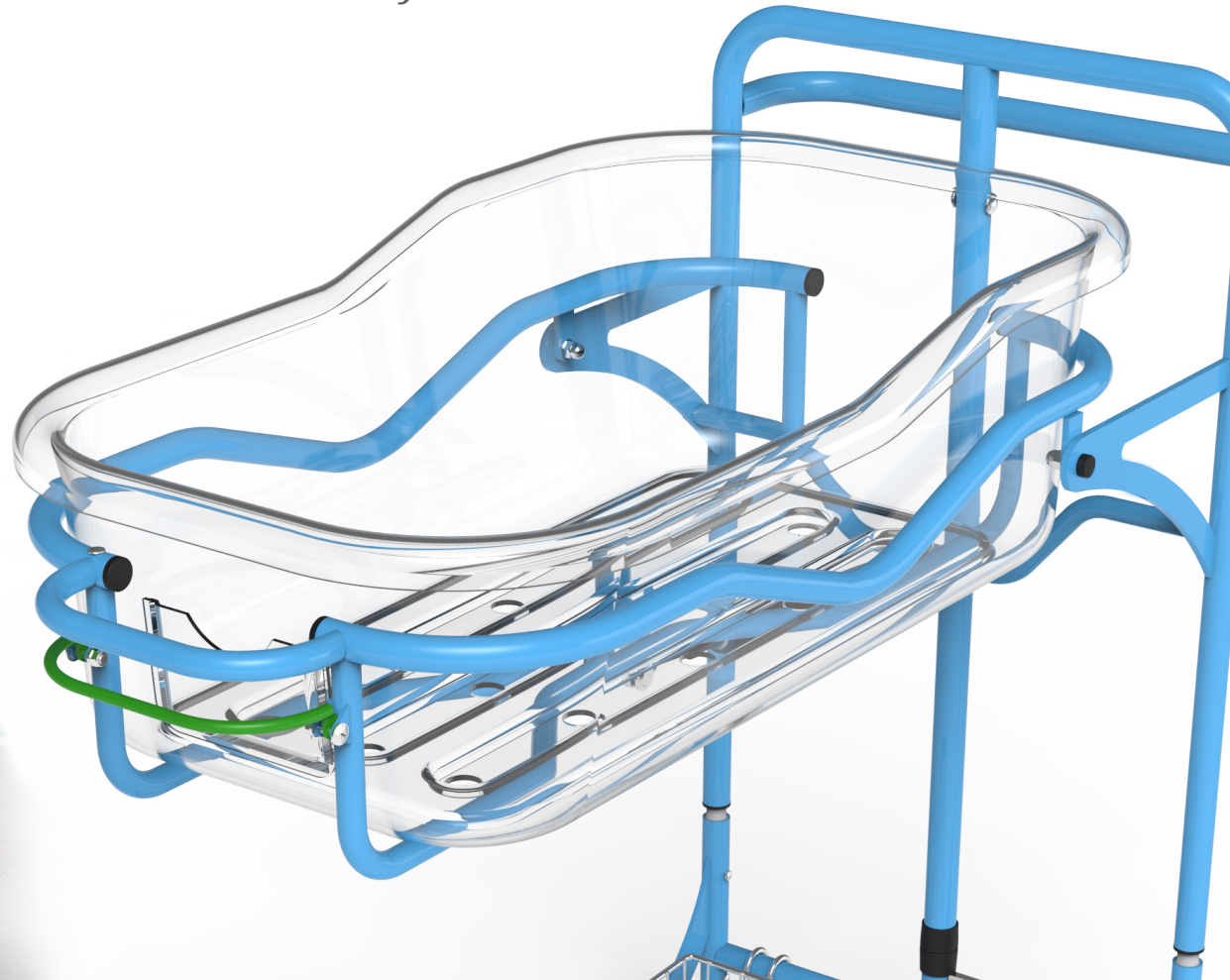

MANIPULATION HANDLE

For easy manipulation and safe transport is newborn bed equipped with manipulation handle.

REMOVABLE TRANSPARENT PLASTIC TUB

Plastic tub is made of PET material which complies with the highest requirements and standards for use in contact with the smallest children. Tub is equipped with holes in lower part. These holes ensure the best air circulation below the mattress.

INTELLIGENT CONSTRUCTION

The bed design offers a wide range of positioning, including height adjustment, which enables the bed to be placed over the mother's bed and set up perfect position for newborn.

TUB LIFTING

You can choose from a wide range of different designs of the newborn beds.

Basic version without tub adjustment or higher versions which enable height adjustment of the tub controlled by hand or foot gas-spring pedal.


WIRE BASKET

Undercarriage is equipped with chromed wire basket for everything you might need for the baby care.

WHEELS

High-quality bearings make transport and handling of the newborn bed much easier and user friendly.

COLOUR SAMPLER

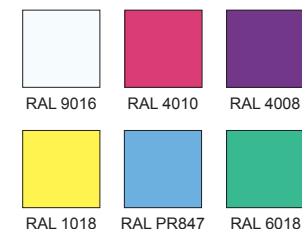

NAMEPLATE

Tub is equipped with transparent pocket where a label with child's name can be simply inserted.

TUB TILTING

The optimum position for the baby may be easily adjusted to several positions even in the basic version by means of a positioning rack. In the higher version the tub may be easily adjusted to any position with one hand by means of a gas spring.



TSE spol. s r. o.
Mánesova 74, 371 52 České Budějovice
T: 386 721 149, F: 386 721 102
medical@tse.cz, www.tse.cz

STANDARD VERSIONS	FIXED HEIGHT		HEIGHT ADJUSTING BY LEG		HEIGHT ADJUSTING BY HAND	
	Version 1	Version 2	Version 3	Version 4	Version 5	Version 6
Tub dimension	75 x 40 cm	75 x 40 cm	75 x 40 cm	75 x 40 cm	75 x 40 cm	75 x 40 cm
Tub height adjustment	✗	✗	✓	✓	✓	✓
Tilt mechanism of tub	mechanical	gas-spring	mechanical	gas-spring	mechanical	gas-spring
Tub tilting (up/down)	✓	✓	✓	✓	✓	✓
Newborn age up to	4 months	4 months	4 months	4 months	4 months	4 months
Maximum safe load	10 kg	10 kg	10 kg	10 kg	10 kg	10 kg
Wire basket	✓	✓	✓	✓	✓	✓
Mattress dimension	74 x 39 cm	74 x 39 cm	74 x 39 cm	74 x 39 cm	74 x 39 cm	74 x 39 cm
Mattress height	4 cm	4 cm	4 cm	4 cm	4 cm	4 cm